

RingJect™

Preloaded Capsular Tension Ring in a **single use injector**

The RingJect™ system is an OPHTEC preloaded Capsular Tension Ring (CTR) in a single use injector. The CTR is prepositioned in the injector and is designed to be self-loading for the surgeon's convenience. The OPHTEC CTR is a high precision medical device for insertion in the capsular bag and made from highly flexible compression molded polymethylmethacrylate (PMMA). Clockwise or counterclockwise insertion is possible. The RingJect™ system is provided in a sterile blister, ready for use.

OPHTEC Capsular Tension Rings are indicated for the stabilization of the capsular bag in the presence of weakened or compromised zonules.

Features & Benefits

Preloaded:	<ul style="list-style-type: none">• Self-loading, single packaging <i>a valuable addition to your surgical armamentarium</i>
Injector:	<ul style="list-style-type: none">• Easy to use• Total control <i>includes action/retraction mechanism</i>• Implantation - clockwise or counter clockwise <i>indicators on injector</i>• Long small tip <i>especially convenient in cases of deep set eyes</i>
Capsular Tension Ring:	<ul style="list-style-type: none">• Circular expansion and stabilization of the capsular bag• Safe IOL centration in eyes with zonular dehiscence• Prevents IOL decentration after capsular shrinking• Stabilized conditions during Phaco-emulsification surgery• Reduced risk of capsular fibrosis• Improves visual acuity when implanted along with premium IOL

RingJect PC375 | CTR Model 275 12/10



Compression:	From 12 mm to 10 mm
Material:	PMMA Flexible Ring
Available in:	Clear PMMA
Overall Ø:	12 mm

RingJect PC376 | CTR Model 276 13/11




Compression:	From 13 mm to 11 mm
Material:	PMMA Flexible Ring
Available in:	Clear PMMA
Overall Ø:	13 mm




Operating instructions

1. Remove the RingJect™ from its blister pack;

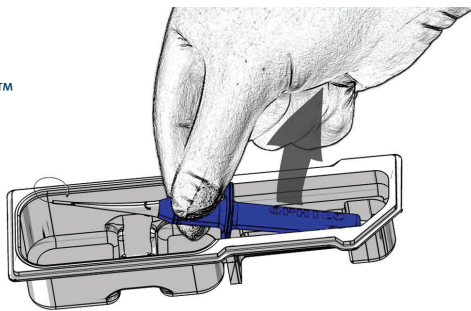


2. Depress the plunger to engage the Capsular Tension Ring, stop pushing after the “click”;

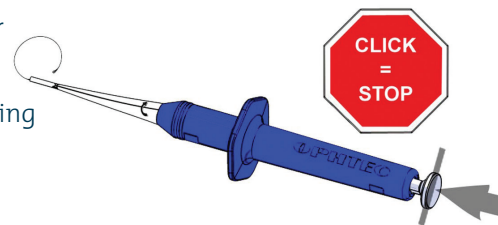


3. Slowly release the plunger, the Capsular Tension Ring will load in the injector;




4. The injector is now loaded and ready for insertion into the eye.

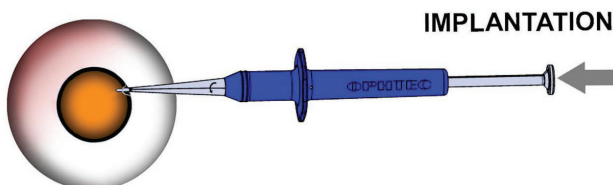

>> Counter clockwise implantation

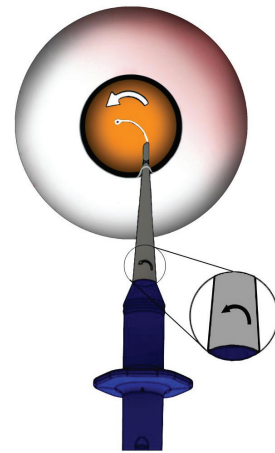

>> Clockwise implantation

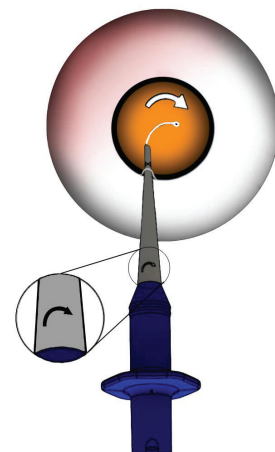

One instrument for clockwise or counter clockwise implantation. The arrow on the base of the tip of the injector shows the direction in which the capsular tension ring will exit the injector.

www.ophtec.com